

Dia Pro®

Across Manuel System







Manuel System

across SYSTEM



AcrossSystem® Card Centrifuge



- 24 Across Gel® Card capacity
- Speed, braking and time optimised for Across Gel® Cards.
- Visual and audible alarm in case of unbalanced loading, power failure or at the end of processing time.
- Remaining processing time and centrifuge speed are shown on the digital display.
- Removable and modifiable impeller and card slots

Ref. No: 900100



AcrossSystem® Card Incubator



- 24 Across Gel® Card capacity
- Two different blocks with 12 Across Gel® Card capacity each working independently.
- Previously adjusted at 37 °C for Across Gel® Cards and time calibrated.
- Remaining processing time and temperature are shown on the digital display designed for two different blocks.
- Visual and audible alarm in case of power failure or at the end of processing time.

Ref. No: 900110

Manuel System



AcrossSystem® Pipette

- Serial pipetting is easier with single intake which allows 8 times 50 ul, 16 times 25 ul or 40 times 10 ul dispensing.
- Easy access to the base of test tubes through pipette tip.
- Single use pipette tips can be disposed quickly and safely.
- Can be applied with both hands thanks to its ergonomic design.

Ref. No: 900120



AcrossSystem® Dispenser

- Fast and reliable usage
- Usage with 100ml or 500ml vials

Ref. No: 900130




AcrossSystem® Workstation

- 12/18/24 Across Gel® Card capacity
- 100ml Across Lisspace
- 500ml Across Lisspace
- 25ml Across Bromelin space
- 8 RBC vial spaces
- 36/54/72 tube areas for preparation of suspension
- Adjustable to work with both hands

Ref. No: 900140

AcrossSystem® Card Reader

Ref. No: 900160

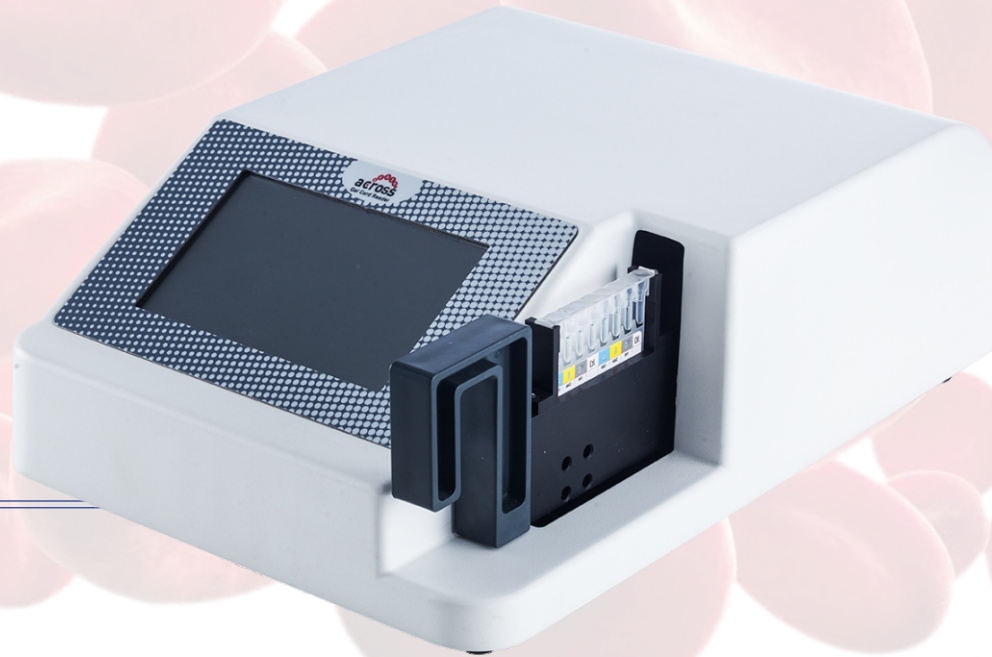

GENERAL FEATURES

- Easy to place in small spaces with a unique, compact design and operating system.
- Across System® Readertype and test type are automatically perceived through the gel card barcodes and the necessary interpretation is presented without user intervention by means of its software.

- Agglutinations in the microtubes are detected through the original display of Across System® Reader by means of its CCD camera. Across System® Test evaluation is performed automatically. It allows user intervention on a touch screen when necessary.
- It is possible to connect Across System® Reader with automation systems dually. Among with test results, Across System® Reader display can also be transferred to automation system when necessary.

Dia Pro®

Tıbbi Ürünler Sanayi ve Ticaret A.Ş.

Merkez / HeadOffice

Barbaros Mah. Tophaneliođlu Cad. Işık Apt No: 70/1 PK 34662 Üsküdar - İstanbul / Türkiye
Tel: +90 (216) 326 70 42 • Faks: +90 (216) 340 16 89 • e-mail: diapro@diapro.com.tr

Fabrika /Factory

İnönü Mah. Gebze Plastikçiler OSB Atatürk Blv. No:25/1 PK 41400 Gebze - Kocaeli / Türkiye
Tel: +90 (262) 751 53 34 • Faks: +90 (262) 751 53 37 • e-mail: diapro@diapro.com.tr



TK1785

www.diapro.com.tr