

Dia Pro®

Across Auto System octo-m







Across Auto System octo-m



Capacity

- 200 Across Gel® Card loading capacity
- 6x16 sample loading capacity
- 12 RBC loading capacity
- 4 spaces for 100ml Across Liss and/or Bromelin
- 3 spaces for 96 well plates
- Incubator with 24 Across Gel® Card capacity
- Centrifuge with 24 Across Gel® Card capacity
- Waste container

TECHNICAL FEATURES

- Enables performing two tests simultaneously via two independent probes.
- Complete walk_away. Across Auto System octo m device performs each step of test procedures as walk-away.
- Across Auto System octo m device provides continuous loading of samples.
- Enables loading various tests simultaneously with its Random Accessmode.
- Both having 24-Card capacity, Incubator and Centrifuge make tests ready in the shortest possible time with its intuitive software.
- An optional RBC cooling block provides longer stability for RBC on the system.
- CE compliant under the IVD Directive of 98/79 EC.



TECHNICAL SUPPORT

Dia Pro or its local distributor provides high quality technical support on national and international market.



Ref. No: 900200


User friendliness of octo-m



- For result evaluation, real card images are processed and stored in archive files onboard the instrument and /or transferred to the hospitals' LIS System. Extensive search functions are available to locate results.
- Prevents mistakes with its dual automation LIS connection.
- By looking at priority samples from a different point of view, Emergency Sample and Stat Sample entry features make a difference.
- Access to result by multi-level password protection. If a second change is made on a test result, the system shows the user ID who made the change together with date and time.
- Provides necessary information on real time loaded tests and their remaining processing times.
- Quick response to required new test panels.
- With the reflex test feature, Weak RH (D) / Lectin A1-Anti H / Antibody titration etc. Tests can be performed automatically.
- Quantitative reports can be accessed according to single or multiple test types within the given dates.
- Remote Access mode
- Multi-language choice. Different operator languages available.
- Average throughput of 360 tests per hour

TEST METHODS

- Forward Grouping
- Forward & Reverse Grouping
- Reverse Grouping
- Rh and Rh Sub-Grouping
- Direct Coombs
- Cross Match
- Antibody Screening
- Antibody Identification
- Antibody Titration
- Weak D / Partial D
- Minor Antigen Tests

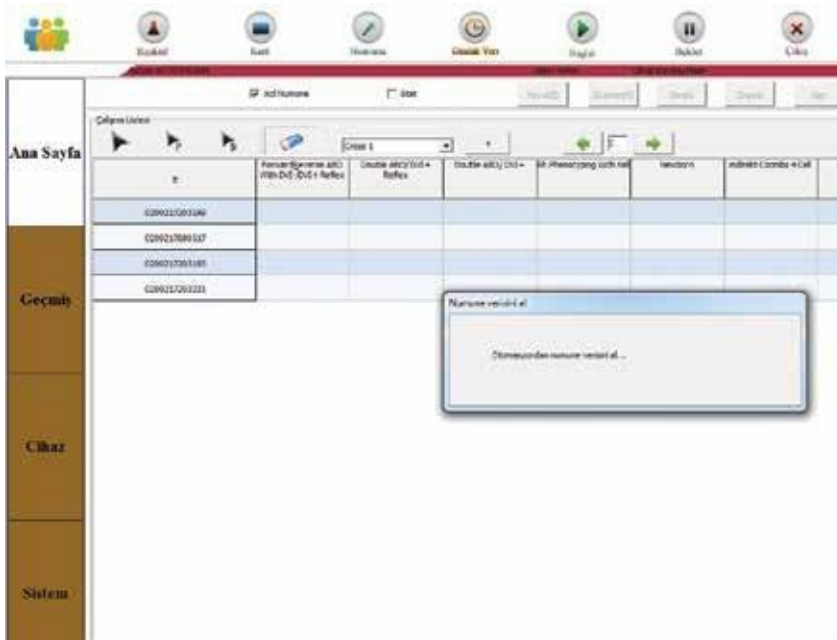

Broad testing parameters
**make easy application
and save time!**



Notlar

Dia Pro®

Tıbbi Ürünler Sanayi ve Ticaret A.Ş.

Merkez / HeadOffice

Barbaros Mah. Tophaneliođlu Cad. Işık Apt No: 70/1 PK 34662 Üsküdar - İstanbul / Türkiye
Tel: +90 (216) 326 70 42 • Faks: +90 (216) 340 16 89 • e-mail: diapro@diapro.com.tr

Fabrika /Factory

İnönü Mah. Gebze Plastikçiler OSB Atatürk Blv. No:25/1 PK 41400 Gebze - Kocaeli / Türkiye
Tel: +90 (262) 751 53 34 • Faks: +90 (262) 751 53 37 • e-mail: diapro@diapro.com.tr



TK1781

www.diapro.com.tr